

β-Catenin clone 14


Instructions For Use

Specification:

β-Catenin is a 92 kDa protein that binds to the cytoplasmic tail of E-Cadherin. The cadherins, transmembrane adhesion molecules, are found with catenins at adherens junctions (zonula adherens). Deletions in the cytoplasmic domain of E-Cadherin which eliminate catenin binding also result in a loss of cell adhesion, indicating that this binding is essential for E-Cadherin function. Although the α- and β-Catenins have been cloned, very little is known about their biochemical roles. However a link between β-Catenin and colon cancer has been described. β-Catenin was found to co-immunoprecipitate with the APC tumor suppressor protein in human colorectal tumor cell lines, as well as in human kidney 293 cells. E-Cadherin, however, was not detectable in these complexes. Thus the APC-Catenin complex may be affecting the transmission of contact inhibition signals and/or the regulation of cell adhesion.

Availability:

Catalog No.	Contents	Volume
ILM 7324 C1	B-Catenin	1,0 ml
ILM 7324 C05	B-Catenin	0,5 ml
ILM 7324 C01	B-Catenin	0,1 ml

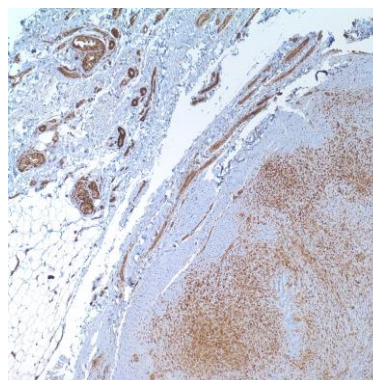

Intended use: For In Vitro Diagnostic Use (IVD)

Reactivity: Human, Chicken, Dog, Mouse, Rat

Clone: 14/Beta-catenin

Species of origin: Mouse

Isotype: IgG1

Control Tissue: Abdomen, fibromatosis of breast, breast carcinoma, prostate, transitional cell carcinoma

Staining: Nuclear, membranous, cytoplasmic

Immunogen: Mouse β-Catenin aa. 571-781

Presentation: 0.05 mg/ml in Aqueous buffered solution containing BSA, glycerol, and ≤0.09% sodium azide.

Application and suggested dilutions:

Pretreatment: Heat induced epitope retrieval in 10 mM citrate buffer , pH6.0, or in 50 mM Tris buffer pH9.5, for 20 minutes is required for IHC staining on formalin-fixed, paraffin embedded tissue sections.

- Immunohistochemical staining of formalin-fixed, paraffin embedded tissue section (dilution up to 1:100-1:200)
- WB (dilution up to 1:250)

The optimal dilution for a specific application should be determined by the investigator.

Note: Dilution of the antibody in 10% normal goat serum followed by a goat anti-mouse secondary antibody based detection is recommended

Storage & Stability: Store at 2-8 °C. Do not use after expiration date printed on the vial.

References:

- 1) Eger A, *J Cell Biol.* 2000; 148(1):173-187
- 2) Lee MSJ *Cell Biol.* 2002; 158(6):1067-1078.
- 3) Ozawa M, *Proc Natl Acad Sci U S A.* 1990; 87(11):4246-4250.
- 4) Persad S, *J Cell Biol.* 2001; 153(6):1161-1173.
- 5) Tateishi K, *J Cell Biol.* 2001;155(4):571-579.