

Alpha Fetoprotein Ab-1 (Clone C3)

Mouse Monoclonal Antibody

Cat. #MS-226-P0, -P1, or -P (0.1ml, 0.5ml, or 1.0ml at 200µg/ml) (Purified Ab with BSA and Azide)

Cat. #MS-226-P1ABX or -PABX (0.1ml or 0.2ml at 1.0mg/ml) (Purified Ab without BSA and Azide)

Description: AFP is normally synthesized in the liver, intestinal tract, and yolk sac of the fetus. Antibody to AFP has been shown to be useful in detecting hepatocellular carcinomas (HCC) and germ cell neoplasms, especially yolk sac tumors.

Mol. Wt. of Antigen: 70kDa

Epitope: Not determined

Species Reactivity: Human, Dog, and Pig. It does not react with cow, cat, mouse, and rat.

Clone Designation: C3

Ig Isotype: IgG_{2a}

Immunogen:

Raised against affinity purified alpha fetoprotein (AFP) from serum of a hepatoma patient

Applications and Suggested Dilutions:

- Immunohistology (Ab-2 is better)
(Ab 1:10-1:20 for 30 min at RT)
- * [No special pretreatment is required for immunohistochemical staining of formalin-fixed, paraffin-embedded tissues]

The optimal dilution for a specific application should be determined by the investigator.

Positive Control:

Fetal liver or hepatocellular carcinoma.

Cellular Localization: Cytoplasmic

Supplied As:

200µg/ml antibody purified from the ascites fluid by Protein A chromatography. Prepared in 10mM PBS, pH 7.4, with 0.2% BSA and 0.09% sodium azide. Also available without BSA and azide at 1mg/ml.

Storage and Stability: Ab with sodium azide is stable for 24 months when stored at 2-8°C. Antibody

WITHOUT sodium azide is stable for 36 months when stored at below 0°C.

Key Reference:

I. Yazova AK; Goussev AI; Poltoranina VS; Yakimenko EF. Human alpha-fetoprotein epitopes as revealed by monoclonal antibodies. Immunology Letters, 1990 Sep, 25(4):325-30.

Limitations and Warranty:

Our products are intended FOR RESEARCH USE ONLY and are not approved for clinical diagnosis, drug use or therapeutic procedures. No products are to be construed as a recommendation for use in violation of any patents. We make no representations, warranties or assurances as to the accuracy or completeness of information provided on our data sheets and website. Our warranty is limited to the actual price paid for the product. NeoMarkers is not liable for any property damage, personal injury, time or effort or economic loss caused by our products.

Material Safety Data:


This product is not licensed or approved for administration to humans or to animals other than the experimental animals. Standard Laboratory Practices should be followed when handling this material. The chemical, physical, and toxicological properties of this material have not been thoroughly investigated. Appropriate measures should be taken to avoid skin and eye contact, inhalation, and ingestion. The material contains 0.09% sodium azide as a preservative. Although the quantity of azide is very small, appropriate care should be taken when handling this material as indicated above. The National Institute of Occupational Safety and Health has issued a bulletin citing the potential explosion hazard due to the reaction of sodium azide with copper, lead, brass, or solder in the plumbing systems. Sodium azide forms hydrazoic acid in acidic conditions and should be discarded in a large volume of running water to avoid deposits forming in metal drainage pipes.

For Research Use Only

Thermo Fisher Scientific
Anatomical Pathology
47777 Warm Springs Blvd.
Fremont, CA 94539, USA
Tel: 1-510-991-2800
Fax: 1-510-991-2826
<http://www.thermo.com/labvision>



Manufactured by:
NeoMarkers
For
Lab Vision Corporation

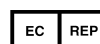

Thermo Fisher Scientific
Anatomical Pathology
93-96 Chadwick Road, Astmoor
Runcorn, Cheshire WA7 1PR, UK
Tel: 44-1928-562600
Fax: 44-1928-562627
Labvision.uk@thermofisher.com